

SMAGEXP: a galaxy tool suite for transcriptomics data meta-analysis

Samuel Blanck^{1,*}, Guillemette Marot^{1,2}

February 26, 2018

¹Univ. Lille Droit et Santé, EA 2694, F-59000 Lille, France

²Inria Lille-Nord Europe, MODAL, F-59000 Lille, France

Abstract

Background: With the proliferation of available microarray and high throughput sequencing experiments in the public domain, the use of meta-analysis methods increases. In these experiments, where the sample size is often limited, meta-analysis offers the possibility to considerably enhance the statistical power and give more accurate results. For those purposes, it combines either effect sizes or results of single studies in a appropriate manner. R packages metaMA and metaRNASeq perform meta-analysis on microarray and NGS data, respectively. They are not interchangeable as they rely on statistical modeling specific to each technology.

Results: SMAGEXP (Statistical Meta-Analysis for Gene EXpression) integrates metaMA and metaRNASeq packages into Galaxy. We aim to propose a unified way to carry out meta-analysis of gene expression data, while taking care of their specificities. We have developed this tool suite to analyse microarray data from Gene Expression Omnibus (GEO) database or custom data from affymetrix[©] microarrays. These data are then combined to carry out meta-analysis using metaMA package. SMAGEXP also offers to combine raw read counts from Next Generation Sequencing (NGS) experiments using DESeq2 and metaRNASeq package. In both cases, key values, independent from the technology type, are reported to judge the quality of the meta-analysis. These tools are available on the Galaxy main tool shed. Source code, help and installation instructions are available on github.

Conclusion: The use of Galaxy offers an easy-to-use gene expression meta-analysis tool suite based on the metaMA and metaRNASeq packages.

Background

Meta-analyses are widely used in medicine and health policy to increase statistical power in studies suffering from small sample sizes. Gene expression experiments are a typical example of such designs. The R packages metaMA and metaRNASeq are dedicated to gene expression microarray and NGS meta-analysis, respectively. While metaMA and metaRNASeq are open source and available on CRAN, they require coding skills in R to perform meta-analysis. Thus, to facilitate the use and the dissemination of these packages, we developed Galaxy wrappers. Galaxy [1, 2, 3] is an open, web-based platform for data intensive biomedical research. It keeps tracks of history and all analyses can be rerun. Galaxy community is very active and a lot of bioinformatics tools are included in Galaxy thanks to a modular system based on XML wrappers. These integrated tools can be shared via the Galaxy toolshed which serves as an appstore.

*to whom correspondence should be addressed

Methods

Overview of R packages integrated into Galaxy

metaMA

Gene expression microarray data meta-analysis can be performed thanks to the metaMA [4] R package. It proposes methods to combine either p-values or moderated effect sizes from different studies to find differentially expressed genes. In our pipeline we only keep the inverse normal method [5] to combine the p-values calculated by limma [6] for each single study.

metaRNAseq

RNA-seq data meta-analysis can be performed thanks to the metaRNASeq [7] R package. It implements two p-value combination techniques : the inverse normal and Fisher methods [8]. Single study p-values are computed with DESeq2 [9].

Differences between metaMA and metaRNASeq

Main differences come from the statistical distributions used to model data and from the manner to treat the genes exhibiting conflicting expression patterns (i.e., under-expression when comparing one condition to another in one study, and over-expression for the same comparison in another study). Usually, microarray data are modelled by Gaussian distributions while NGS data are modelled by Negative Binomial distributions. As explained in [4] and [7], the trick which consists to use one-tailed p-values for each single study before combination in metaMA avoids directional conflicts. In metaRNASeq, this trick can not be used, which necessitates a post-hoc identification of conflicts, step which is also proposed in metaRNASeq.

Description of Galaxy tools

SMAGEXP tool suite offers two distinct gene expression meta-analysis functionalities : one dedicated to microarray data meta-analysis and one dedicated to RNAseq data meta-analysis (see Table 1 and figure 1).

Table 1: Summary of tools inputs and outputs.

Tool	Input	Output
GEOQuery	GEO database ID	Rdata object and .cond file
QCNormalization	Raw .CEL affymetrix files and .conf file	Rdata object and plots
Import custom data	Expression data in tabular .txt format	Rdata object and plots
Limma analysis	Rdata object from GEOQuery or QCNormalization	Rdata Object and HTML report
Microarray data meta-analysis	Rdata objects from Limma analyse	HTML report
RNA-seq data meta-analysis	Results text files from galaxy deseq2 tool	HTML report

Microarray data meta-analysis

GEOQuery tool GEOQuery tool fetches microarray data directly from GEO database, based on the GEOQuery [10] R package. Given a GSE accession ID, it returns an Rdata object containing the data and a text file (.cond file) summarizing the conditions of the experiment. The .cond file is a text file containing one line per sample in the experiment. Each line is made of 3 columns:

- Sample ID

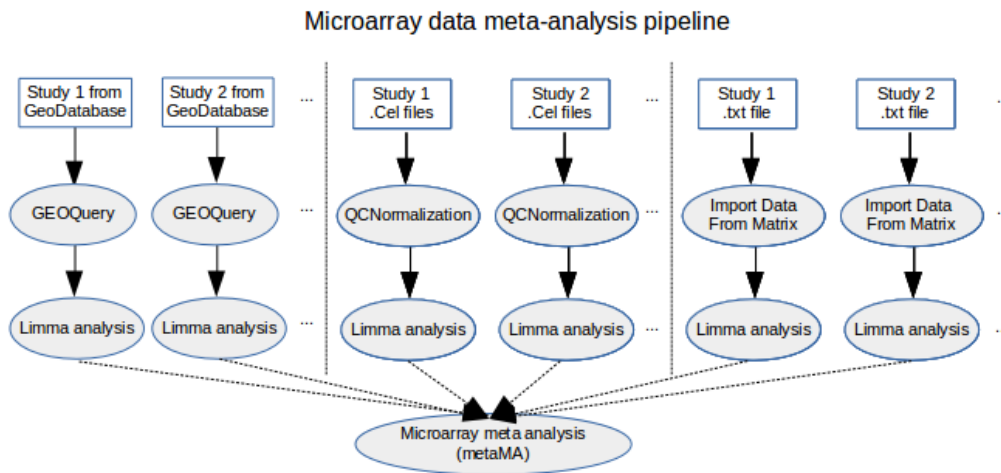


Figure 1: Overview of the tools from microarray data meta-analysis pipeline integrated within Galaxy.

- Condition of the biological sample
- Description of the biological sample

```

GSM342582.CEL    tumor    GSM342582_Tongue_040
GSM342583.CEL    normal   GSM342583_Tongue_041
GSM342584.CEL    tumor    GSM342584_Tongue_041
GSM342585.CEL    normal   GSM342585_Tongue_042
GSM342586.CEL    tumor    GSM342586_Tongue_042
GSM342587.CEL    normal   GSM342587_Tongue_043
  
```

Figure 2: Example of .cond file.

Column names are optional and only the columns order matters. As the GEO dataset should already have been normalized, the GEOQuery tool does not perform any normalization method, apart from an optional log2 transformation.

QCNormalization tool It is possible to analyze .CEL files from affymetrix gene expression microarray. The QCNormalization tool offers to ensure the quality of the data and to normalize them. Several normalization methods are available :

- rma normalization
- quantile normalization + log2
- background correction + log2
- log2 only

This tool generates several quality figures : microarray images, boxplots and MA plots. It also outputs an Rdata object containing the normalized data for further analysis with the limma analysis tool.

Import custom data tool This tool imports data stored in a tabular text file. Column titles (chip IDs) must match the IDs of the .cond file. A few normalization methods are proposed, but it is possible to skip the normalization step, by choosing "none" in the normalization methods options. Therefore this tool is of special interest when the input dataset has been previously normalized.

This tool also generates boxplots and MA plots and outputs an Rdata object containing the data for further analysis with the limma analysis tool.

Limma analysis tool The Limma analysis tool performs single analysis either of data previously retrieved from GEO database or normalized affymetrix .CEL files data. Given a .cond file, it runs a standard limma differential expression analysis. The user choose two conditions extracted from the .cond file (see Figure 3). It generates boxplots for rough quality control of normalization, p-value histograms to ensure that statistical hypotheses are not violated and a volcano plot to quickly identify the most-meaningful changes. This tool also outputs a table summarizing the differentially expressed genes and their annotations. Genes are sorted by ascending Benjamini-Hochberg adjusted p-value, and annotations are retrieved via GEO database. This list of genes can be exported to excel or to csv format. This table is sortable and requestable. Furthermore it is possible to expand each row to display extended annotations informations, including hypertext links to the National Center for Biotechnology Information (NCBI) gene database. Finally, this tool outputs an Rdata object to perform further meta-analysis and a tabular file containing the all results and annotations of the differential analysis.

The screenshot shows the Galaxy interface for the Limma analysis tool. The title bar reads "Limma analysis Performs gene expression analysis thanks to limma (Galaxy Version 1.0.0)". Below the title bar, there are several input fields and a button:

- RData:** A dropdown menu showing "42: GEOQuery RData of GSE3524".
- Conditions:** A dropdown menu showing "41: conditions of GSE3524".
- condition 1:** A dropdown menu showing "series of 16 tumors".
- condition 2:** A dropdown menu showing "series of 4 normals".
- number of top genes:** A text input field containing "1000".
- Execute:** A blue button with a checkmark icon and the text "Execute".

Figure 3: limma analysis tool form.

Microarray data meta-analysis tool The meta-analysis relies on the metaMA R package. Prior to the meta-analysis itself, a pre-processing is made in order to ensure compatibility between several sources of data. In fact, data could come from different types of microarrays. First, we list the Entrez gene ID corresponding to each probe of each dataset. Next, we keep the probes corresponding to the genes which are shared by all the experiments of the meta-analysis. Then, for each dataset, we merge the microarray probes originating from the same Entrez gene ID by computing their mean. Note that the merging of different technologies induces a loss of

information and might generate several conflicts as probes do not necessary reflect the same biological reality. Finally, the p-value combination method of metaMA is run on the merged dataset. It generates a Venn Diagram summarizing the results of the meta-analysis, and a list of indicators to evaluate the quality of the performance of the meta-analysis :

- DE : Number of differentially expressed genes
- IDD (Integration Driven discoveries) : number of genes that are declared differentially expressed in the meta-analysis that were not identified in any of the single studies alone
- Loss : Number of genes that are identified differentially expressed in single studies but not in meta-analysis
- IDR (Integration-driven Discovery Rate) : corresponding proportion of IDD
- IRR (Integration-driven Revision) : corresponding proportion of Loss

It also outputs a fully sortable and requestable table, with gene annotations and hypertext links to NCBI gene database.

RNA-seq data meta-analysis

The RNA-seq data meta-analysis tool relies on the deseq2 galaxy tool analysis results. Given several text file resulting from the deseq2 tool, the metaRNAseq tool performs a meta-analysis, generates the list of differentially expressed genes, and outputs the DE, IDD, Loss, IDR and IRR indicators.

Application

Microarray meta-analysis example

SMAGEXP was applied to two GEO datasets identified with the following IDs : GSE3524 and GSE13601. These two datasets contain human oral squamous cell carcinoma (SCC) data. See Figure 4 for an overview of the workflow of this analysis.


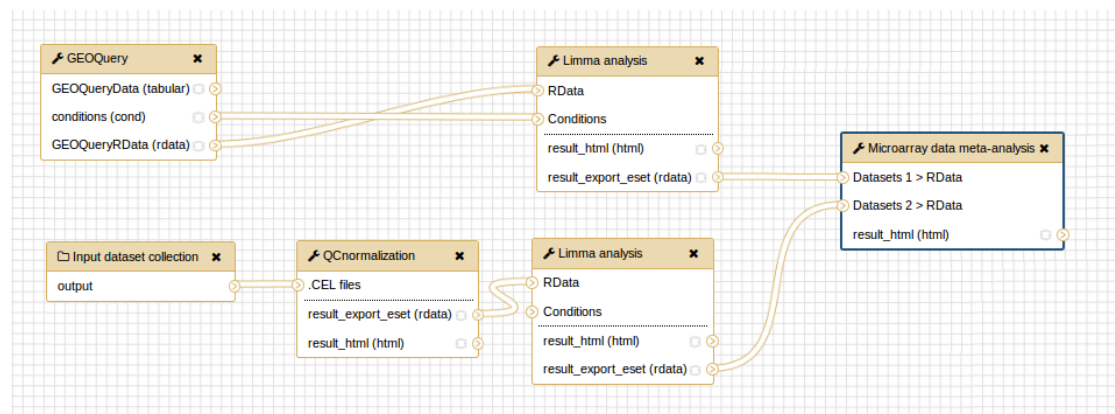


Figure 4: Exemple of a galaxy workflow for microarray meta-analysis.

First, we fetch data from the GSE3524 using the GEOQuery tool (with parameter "log2 transformation" = auto). Then we launch the limma analysis, using the output from the GEOquery tool. It generates an Rdata output, which will be usefull for the meta-analysis. Results can be seen on Figure 5 and Figure 6

Boxplots

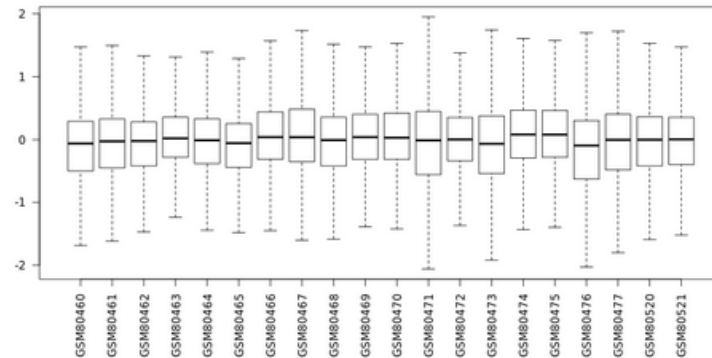

P-value histogram and Volcano plot

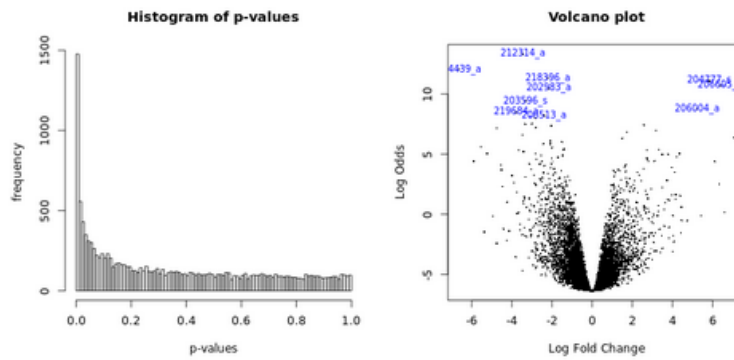


Figure 5: limma analysis tool output plots.

ID	adj_P_Val	P_Value	t	B	logFC	Gene_symbol	Gene_title	Gene_ID	Chromosome_annota...	GO_Function_ID
212314_at	3.3e-06	2.3e-10	-11.19	13.3288	-3.46	SEL1L3	sel-1 suppressor of lin...	23231	Chromosome 4, NC_00...	
204439_at	7.1e-06	1.0e-09	-10.31	12.0465	-6.64	IFI44L	interferon-induced pr...	10964	Chromosome 1, NC_00...	
218396_at	1.0e-05	2.3e-09	-9.85	11.3259	-2.20	VPS13C	vacuolar protein sorti...	54832	Chromosome 15, NC_0...	
204777_s_at	1.0e-05	2.9e-09	9.71	11.1038	5.81	MAL	mal, T-cell differentiat...	4118	Chromosome 2, NC_00...	GO:0015267///GO:000...
Gene Symbol: MAL										
Gene Title: mal, Tcell differentiation protein										
GO Function ID: GO:0015267 , GO:0008289 , GO:0016505 , GO:0005515 , GO:0019911										
206605_at	1.3e-05	4.5e-09	9.47	10.7170	6.39	ENDOU	endonuclease, polyU-s...	8909	Chromosome 12, NC_0...	GO:0003723///GO:000...
202983_at	1.4e-05	6.0e-09	-9.32	10.4674	-2.17	HITF	helicase-like transcrip...	6596	Chromosome 3, NC_00...	GO:0005524///GO:001...
203596_s_at	3.9e-05	1.9e-08	-8.70	9.4094	-3.31	IFIT5	interferon-induced pr...	24138	Chromosome 10, NC_0...	GO:0003723///GO:004...
206004_at	7.0e-05	4.0e-08	8.33	8.7509	5.21	TGM3	transglutaminase 3	7053	Chromosome 20, NC_0...	GO:0005509///GO:000...
219684_at	8.3e-05	5.3e-08	-8.18	8.4819	-3.80	RTP4	receptor (chemosenso...	64108	Chromosome 3, NC_00...	GO:0005515
206513_at	1.0e-04	7.1e-08	-8.04	8.2131	-2.39	AIM2	absent in melanoma 2	9447	Chromosome 1, NC_00...	GO:0003690///GO:004...

Show 10 entries
Showing 1 to 10 of 1,000 entries

Figure 6: limma analysis tool : table of top 10 genes for GSE3524 dataset.


Secondly, the same kind of analysis is run from raw .CEL files. We choose to keep six .CEL files from the GSE13601 dataset (IDs from GSM342582 to GSM342587). Quality control and normalization is done thanks to the QCnormalization tool. Then, as previously, the limma

analysis tool is run to generate a HTML report and an Rdata output.

Run a metaMA analysis

To run the microarray meta-analysis tool, we only need the Rdata output of each single study, generated by the limma analysis tool. It generates a Venn diagram to compare the results of each study with the meta-analysis. It also outputs several indicators as described in the description of the tool (see Figure 7). As for the limma tool, annotated expressed genes are displayed in a table which can be ordered and requested.

Venn diagram

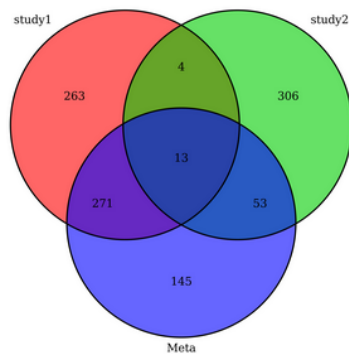

Summary

DE	IDD	Loss	IDR	IRR
482	145	573	30.08	62.97

DE : Number of differentially expressed genes
 IDD (Integration Driven discoveries) : number of genes that are declared DE in the meta-analysis that were not identified in any of the individual studies alone
 Loss : Number of genes that are identified DE in individual studies but not in meta-analysis
 IDR (Integration-driven Discovery Rate) : corresponding proportions of IDD
 IRR (Integration-driven Revision) : corresponding proportions of Loss

Figure 7: Venn diagram and summary of microarray data meta-analysis tool results

RNA-seq data meta-analysis example

The RNA-seq data meta-analysis tool relies on deseq2 results (see Figure 8).

It outputs a Venn diagram and the same indicators as in the microarray data analysis tool for both Fisher and inverse normal p-values combinations. It also generates a text file containing summarization of the results of each single analysis and meta-analysis. Potential conflicts between single analyses are indicated by zero values in the "signFC" column.

1	2	3	4	5	6	7
tag4997	3120.17146561277	3.77231311375556	0.171247901593321	22.0283756977883	1.5399715529456e-107	2.82938973422694e-103
tag9661	869.701524875407	4.15091605808289	0.198823712087783	20.8773692760057	8.59982353619356e-97	7.90022789152422e-93
tag2250	7662.23983112397	-2.90777142681232	0.14982162124608	-19.4082229429112	6.57624117645105e-84	4.02750930449784e-80
tag345	4959.42596250046	2.44621038150424	0.130388652619725	18.7609146375533	1.57663046421531e-78	7.24185787975697e-75
tag9488	6632.50141736004	-2.61596039303527	0.141457108109693	-18.4929582400817	2.35285730311673e-76	8.64580944603275e-73
tag4264	3215.03391406763	-3.05128926499533	0.165825317089365	-18.4006237319658	1.29861356542041e-75	3.51102594645424e-72

Figure 8: Header of a deseq2 results file

1	2	3	4	5	6	7	8
"ID"	"DE.DE.deseq_results1.txt"	"DE.DE.deseq_results2.txt"	"DE.DE.fishercomb"	"DE.DE.invnorm"	"FC.deseq_results1.txt"	"FC.deseq_results2.txt"	"signFC"
"tag1"	0	0	0	0	-0.0514474564769571	-0.0996641983598344	-1
"tag10"	0	0	0	0	-0.341968584114381	0.0306325127845984	0
"tag100"	0	0	0	0	0.218292477256233	-0.0691816856430162	0
"tag1000"	0	1	1	1	0.363870914458309	0.346232782219452	1
"tag10000"	0	0	0	0	-0.0805354648414234	-0.126612590668799	-1
"tag10001"	0	0	0	0	0.0717225595372206	0.120161813963893	1
"tag10002"	0	0	0	0	-0.114847586098109	0.0123268534412585	0
"tag10003"	0	0	0	0	0.0217198765889001	-0.0312857938137718	0
"tag10004"	0	0	0	0	-0.173677437216902	-0.09142089339658	-1
"tag10005"	0	0	0	0	-0.04625232644183	-0.0193066755866328	-1
"tag10006"	0	0	0	0	-0.16000054078948	-0.0840294513941694	-1
"tag10007"	0	0	0	0	0.612687034527193	-0.0642077097755686	0
"tag10008"	NA	NA	0	0	0.0792311974558976	-0.000228569599220305	0
"tag10009"	1	1	1	1	-0.962744832164662	1.13707735403279	0

Figure 9: Header of a metaRNAseq results file

Conclusion

We developed SMAGEXP, a toolsuite dedicated to gene-expression data meta-analysis. This toolsuite proposes quality controls, single analyses and meta-analyses of microarray and RNA-seq data, suggesting appropriate pipelines for each type of data. It delivers fully annotated results of differentially expressed genes, exportable in several usual formats. Integrated into Galaxy, SMAGEXP is easy to use for biologists and life scientists. R packages metaMA and metaRNAseq thus inherit reproducibility and accessibility support from Galaxy. SMAGEXP is available on the Galaxy main toolshed [11].

Source code is available on github at : <https://github.com/sblanck/smagexp>.

Furthermore, thanks to Docker, we made these Galaxy tools and their dependencies easy to deploy. A fully dockerized instance of Galaxy containing SMAGEXP is available at : <https://hub.docker.com/r/sblanck/galaxy-SMAGEXP/>

References

- [1] J. Goecks, A. Nekrutenko, J. Taylor, and T. G. Team, "Galaxy: a comprehensive approach for supporting accessible, reproducible, and transparent computational research in the life sciences," *Genome Biol*, vol. 11, no. 8, p. R86, 2010.
- [2] D. Blankenberg, G. V. Kuster, N. Coraor, G. Ananda, R. Lazarus, M. Mangan, A. Nekrutenko, and J. Taylor, "Galaxy: A web-based genome analysis tool for experimentalists," *Current protocols in molecular biology*, pp. 19–10, 2010.

- [3] B. Giardine, C. Riemer, R. C. Hardison, R. Burhans, L. Elnitski, P. Shah, Y. Zhang, D. Blankenberg, I. Albert, J. Taylor, W. C. Miller, W. J. Kent, and A. Nekrutenko, “Galaxy: a platform for interactive large-scale genome analysis,” *Genome research*, vol. 15, no. 10, pp. 1451–1455, 2005.
- [4] G. Marot, J.-L. Foulley, C.-D. Mayer, and F. Jaffrezic, “Moderated effect size and p-value combinations for microarray meta-analyses,” *Bioinformatics*, vol. 25, no. 20, pp. 2692–2699, 2009.
- [5] L. Hedges and I. Olkin, *Statistical Methods for Meta-Analysis*. London: Academic Press, 1985.
- [6] M. E. Ritchie, B. Phipson, D. Wu, Y. Hu, C. W. Law, W. Shi, and G. K. Smyth, “limma powers differential expression analyses for RNA-sequencing and microarray studies,” *Nucleic Acids Research*, vol. 43, no. 7, p. e47, 2015.
- [7] A. Rau, G. Marot, and F. Jaffrézic, “Differential meta-analysis of rna-seq data from multiple studies,” *BMC Bioinformatics*, vol. 15, no. 1, pp. 1–10, 2014.
- [8] R. A. Fisher, *Statistical Methods for Research Workers*. Edinburgh: Oliver and Boyd, 1932.
- [9] M. I. Love, W. Huber, and S. Anders, “Moderated estimation of fold change and dispersion for rna-seq data with deseq2,” *Genome Biology*, vol. 15, p. 550, 2014.
- [10] S. Davis and P. Meltzer, “Geoquery: a bridge between the gene expression omnibus (geo) and bioconductor,” *Bioinformatics*, vol. 14, pp. 1846–1847, 2007.
- [11] D. Blankenberg, G. Von Kuster, E. Bouvier, D. Baker, E. Afgan, N. Stoler, J. Taylor, and A. Nekrutenko, “Dissemination of scientific software with galaxy toolshed,” *Genome Biology*, vol. 15, no. 2, pp. 1–3, 2014.