

Homemade kit for demonstrating Barkhausen Effect

Shantanu Shakya¹ and Navinder Singh²

¹*Department of Physics, Indian Institute of Space Science and Technology(IIST), Thiruvananthapuram, Kerala - 695547*

²*Theoretical Physics Division, Physical Research Laboratory, Navrangpura, Ahmedabad - 380009*

This paper presents an innovative and cost-effective approach to understanding the Barkhausen effect through the design and implementation of an educational kit. The Barkhausen effect, characterized by Barkhausen noise (BN) during magnetization changes in soft magnetic materials, is explored for its application in probing hysteresis properties and magnetization dynamics. The study investigates scaling properties, categorizing ferromagnetic materials based on scaling exponents. The primary contribution is the introduction of a practical and accessible kit for hands-on Barkhausen Effect demonstrations, revolutionizing the educational experience. This kit enables students to not only comprehend the intricacies of BN but also calculate the scaling constant (τ) for Soft Iron samples. The paper demonstrates the successful construction of the kit, its signal amplification capabilities, and data collection accuracy, showcasing its potential for widespread educational use.

I. INTRODUCTION

Barkhausen in 1919 noticed that changes in the magnetization of an iron sample induced a special form of time dependent electric potential in a pickup coil wound over the sample. Upon amplification, this led to the generation of audible "noise" as depicted in Figure 1. This phenomenon is widely known as "Barkhausen noise" [1]. Barkhausen suggested its association with irreversible alterations in magnetization. When exposed

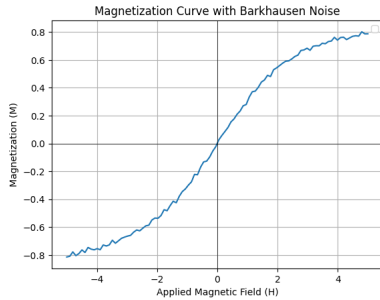


FIG. 1. Barkhausen Noise in magnetization curve

to a magnetic field undergoing a gradual temporal shift, such as the slow movement of a manually-driven magnetic yoke, a usual magnetic material like iron produces irregular signal in a pick-up coil. This irregular acoustic manifestation contrasts sharply with the systematic pattern of the applied magnetic field and is detectable through a microphone. This marked the initial indirect confirmation of the presence of magnetic domains, a concept theorized a few years earlier by Weiss [2].

Soon after its discovery, it became apparent that Barkhausen noise (BN) could serve as a valuable tool for delving into and understanding the magnetization dynamics inherent in soft magnetic materials, provid-

ing insights into specific hysteresis properties. While the measurement of BN is relatively straightforward, its interpretation presents a greater challenge. This complexity stems from the stochastic nature of domain wall (DW) motion, which progresses through sudden jumps or avalanches during slow magnetization. These dynamics are notably influenced by material microstructure, along with additional factors such as the demagnetizing field, external stress, and others [3].

II. EXPERIMENTAL IDEA

In this section, we will delve into the details of the experiment, including the instruments employed, and the associated circuitry. We will explore both the auditory aspect of detecting Barkhausen noise (BN) and the methods for its measurement.

BN is generated by alterations in the magnetization of a ferromagnetic material. This change in magnetization is often attributed to the presence of impurities and defects within the material, which can act as obstacles, pinning the movement of domain walls (as shown in figure 2). When we apply an external magnetic field and increase it beyond a critical value denoted as H_c , a pivotal event known as the "depinning transition" occurs. During this transition, the pinned domain walls are suddenly released or depinned, resulting in a significant jump in the material's magnetization curve. This jump can be detected using a pickup coil that operates on the principles of "Faraday's law of electromagnetic induction".

According to Faraday's law, any alteration in magnetic flux within the coil's vicinity induces a corresponding voltage at the coil's terminals. This induced voltage

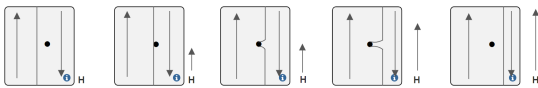


FIG. 2. Pinning of domain wall due to impurities

is the key to our ability to both hear and measure the Barkhausen noise in our experiment.

$$V = N * A * \frac{dB}{dT} \quad (1)$$

V- voltage induced in the coil

N- number of the turns of the pick-up coil

A- cross section area of the coil

B- magnetic field

When deliberately and gradually adjusting the magnetic field strength, either increasing or decreasing, the coil induces a prominent voltage signal characterized by a large amplitude and low frequency. In contrast, the rapid Barkhausen events, characterized by much smaller amplitudes, contribute to a signal with higher frequency. This high-frequency signal is superimposed upon the slower and larger signal. Our aim is to filter out the BN noise voltage and reject the slower and larger signals.

Figure 1 of reference 4[8] illustrates the magnetization curve of a NiFe film [8]. Initially appearing smooth when recorded under the influence of a slowly varying external magnetic field, a closer inspection with magnification reveals non-uniform changes in magnetization. This non-uniformity is marked by the presence of high-frequency, low-amplitude Barkhausen noise, which is superimposed on the main signal.

To extract the Barkhausen noise from the recorded data, we adopt a two-step approach. Initially, a high-gain amplifier is employed to boost the amplitude of the signal. Following this amplification, we implement suitable filtering techniques to isolate the Barkhausen noise. This is accomplished using band-pass or high-pass filters, effectively eliminating lower-frequency components associated with slower and larger signals (refer to Figure 3). Frequencies below 300 Hz are filtered out.[8].

This approach enables us to discern and analyze the Barkhausen noise, a phenomenon of particular interest for further investigation and understanding.

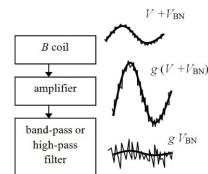


FIG. 3. This is how we are extracting the signal[9]

III. EXPERIMENTAL SET-UP AND PROPOSAL OF KIT

To conduct this experiment, we devised a straightforward homemade kit, from which we extracted the Barkhausen Noise. The kit comprises the following components:

- Ferromagnetic Sample
- Pick-up coil
- Speaker
- Permanent Magnet
- Audio Amplifier

A. Ferromagnetic Sample

In this experiment, we are utilizing a cylindrical soft iron core with a diameter of 19.5 mm and a length of 71 mm. Soft iron is an easily accessible ferromagnetic material characterized by high permeability and low coercivity, meaning it does not retain its magnetization when the magnetic field is removed. This property allows for easy magnetization and demagnetization, and it exhibits minimal hysteresis loss, making it the ideal material for our experiment, where we require repeated switching of the magnetic field to accurately observe Barkhausen noise (BN).

In soft Iron core of not very high purity which is used in standard applications such as in concrete structures etc. These impurities are the primary cause of Barkhausen Noise. Non-magnetic impurities, such as carbon, sulfur, phosphorous, and silicon, are present, along with magnetic impurities like nickel, cobalt, and manganese. Additionally, the material contains grain boundary defects, dislocations within its crystal structure (which can impede the movement of domain walls), as well as inclusions and surface defects.

B. Pick-up Coil

The pickup coil serves the crucial function of detecting and converting alterations in magnetic flux into electrical signals. It plays a pivotal role in various devices where the conversion of magnetic fields into electrical signals is essential for purposes such as detection, measurement, or transmission. The operation of a pickup coil is based on Faraday's law of electromagnetic induction. The magnitude of the voltage signal it produces depends on several factors, including the strength of the magnetic field, the number of turns within the pickup coil, and the rate at which the magnetic field changes.

For our experiment, we have specifically designed a four-layered pickup coil using enameled copper wire with a total of 300 turns. The copper wire employed in this coil has a cross-sectional diameter of 0.5 mm. This coil is wound around the sample (Soft Iron core), enabling it to effectively capture the Barkhausen noise that we are studying in our experiment.



FIG. 4. Pickup coil wound over soft iron (sample)

C. Speaker

A speaker is an electroacoustic transducer designed to transform electrical signals into sound waves. It functions through the interaction between a permanent magnet and an electromagnet that are affixed to a diaphragm. This interaction results in the diaphragm's movement, which corresponds to the incoming electrical signal, ultimately producing sound.

In our experimental setup, crafted for demonstra-

tion purposes, we have incorporated an 8 Ω impedance speaker with a power rating of 0.25 W. This speaker enables us to audibly perceive the faint yet distinct Barkhausen noise (BN) when one listens closely to it.



FIG. 5. speaker used for this experiment

D. Magnet

Certainly, for this demonstration, an external magnetic field is essential for the sample. We are utilizing a "ring-shaped" permanent ferrite magnet, which we retrieved from old speakers. This ferrite magnet can generate a magnetic field with a range spanning from a few milli-tesla (mT) to sub tesla (T), depending on orientation and distance from the magnet.

Ferrite magnets are the most commonly used type of magnets in speakers. They are typically composed of a mixture of iron oxide and other materials, often strontium or barium. These magnets are not only cost-effective but also possess favorable magnetic properties.

When we bring the magnet in proximity to the sample, the sample becomes magnetized, and this magnetization process is accompanied by the generation of Barkhausen noise, which we can audibly perceive.

We are using two ferrite ring magnets for our demonstration. The larger one has a diameter of 110 mm and a height of 15mm, while the smaller one has a diameter of 80 mm and a height of 14mm. These magnets are combined and positioned as illustrated in the figure 6 for our experiment.



FIG. 6. Magnets used.

E. Audio Amplifier

An audio amplifier is an electronic device designed to increase the amplitude of low-power electronic audio signals within the frequency range of approximately 20 to 20,000 Hz. These signals can originate from various sources, such as a radio receiver or an electric guitar pickup, and the amplifier's purpose is to elevate them to a level suitable for driving speakers or headphones.

To amplify the Barkhausen noise (BN) in our experiment, we have developed a cost-effective "Audio Amplifier" utilizing the LM386 integrated circuit (IC). This IC offers adjustable gain and boasts low power consumption. We are able to power this amplifier efficiently using a 9V battery.

1. LM386

The **LM386** is a 8 pin IC, used for designing power amplifier which consume very low power, have adjustable gain of any value from 20 to 200, which could be done by connecting appropriate capacitor between pin 1 and 8 (see Figure 7) [10].

*For more information, see reference on IC LM386[10]

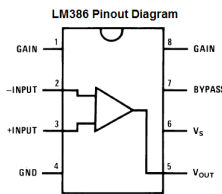


FIG. 7. LM386. For details refer to: www.learningaboutelectronics.com/Articles/How-to-connect-a-LM386-audio-amplifier-chip

F. Amplifier Circuit

Figure 8 illustrates the amplifier circuit constructed using the LM386 IC. This circuit is designed to enhance the Barkhausen noise received from the pickup coil, with pin 3 serving as the input. To optimize the amplification process, a 220Ω resistor is connected in series with pin 3. Additionally, a $10\mu\text{F}$ capacitor is positioned between pins 1 and 8 to achieve maximum gain. A "bypass" capacitor with a value of $100\mu\text{F}$ is connected through pin 7, and another $100\mu\text{F}$ capacitor is linked to the power supply to ground any AC components within the power supply.

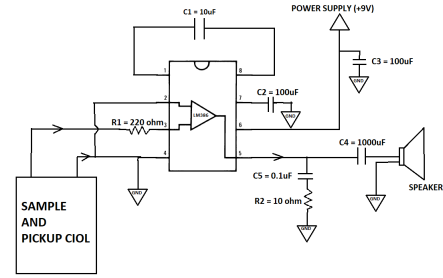


FIG. 8. Circuit designed for this kit

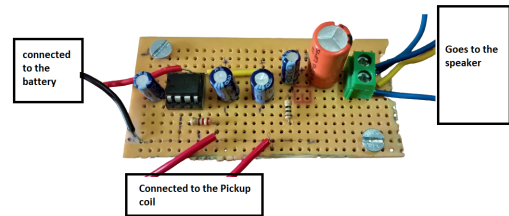


FIG. 9. Circuit made on general purpose PCB

At the output, which is pin 5, a $1000\mu\text{F}$ capacitor is utilized to block the DC signal (by the motion of hands) from reaching the speaker. Furthermore, a low-pass filter is incorporated in this section to eliminate any high-frequency noise (specially MHz noise), particularly in the megahertz range, that may be present along with the BN signal.

Components are mounted on a piece of plywood, and made a very simple kit. (see figure 10)

It turns out that the total cost to make it is only around Rs 300/-.

IV. RESULTS

Initially, we successfully constructed a cost-effective kit (refer to Figure 10), utilizing components detailed in Section III. This kit, as demonstrated in a relevant video accessible as supplementary material at <https://tinyurl.com/BNdata>, effectively captures the Barkhausen Noise, producing clear auditory signals through the integrated speaker. Subsequently, we directed the resulting signal to the *Keithley 2450 sourcemeter*, a versatile instrument serving as a digital multimeter equipped with built-in ammeter, voltmeter, and ohmmeter functionalities. Notably, this instrument boasts a remarkable resolution of 10 nV for the voltmeter

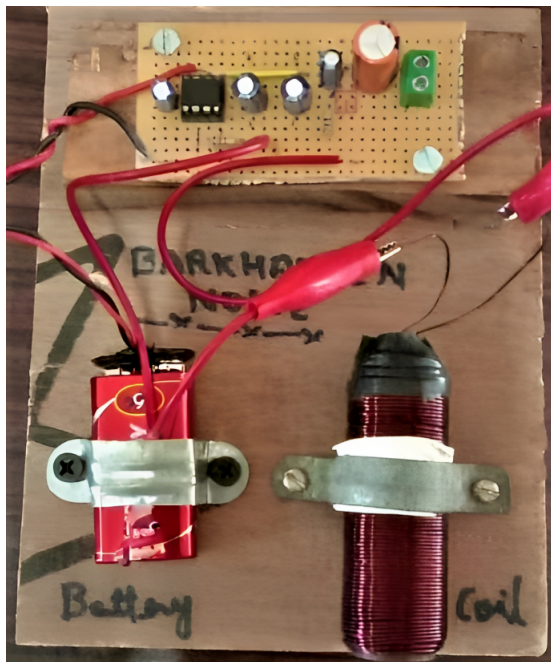


FIG. 10. Barkhausen Experiment kit designed by us.

and 10 fA for the ammeter, featuring 6.5 significant digits and a reading capability ranging from 1 reading per second to 1700 readings per second.

The recorded waveforms are visually presented in Figure 11. Upon bringing the magnet close to the sample,

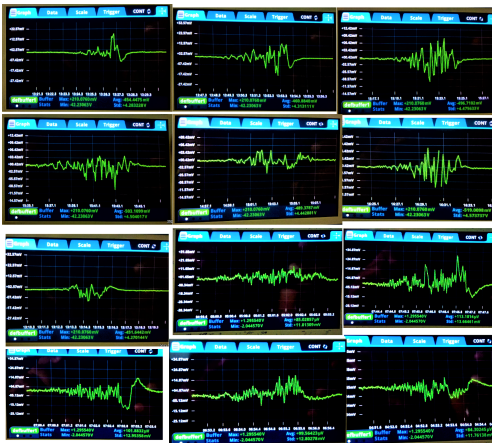


FIG. 11. Barkhausen Avalanche waveforms

the resulting signals are observed in millivolt range, as depicted in the above figure. This indicates the effective performance of our kit, adept at amplifying the microvolt Barkhausen signal into the millivolt range. Subsequently, we conducted measurements to determine the jump size (s) and the frequency ($P(s)$) of the avalanches. The tail part of the measured data was plotted and fit-

ted with a power-law model (refer to Figure 12) given by $P(s) = as^{-\tau} + b$, where 'a' and 'b' represent scaling constants. The obtained τ value is 1.303, with $a=799.139$ and $b=-4.834$.

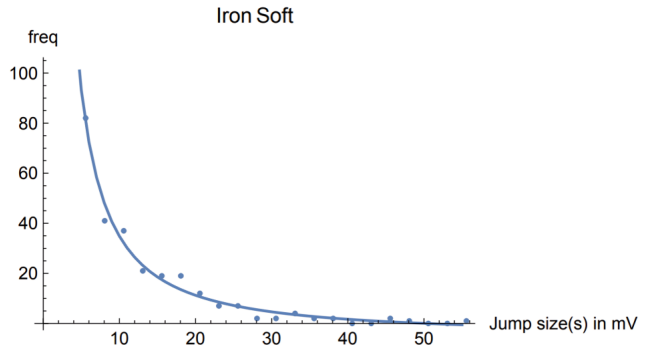


FIG. 12. Curve fitting of observed Barkhausen data.

V. STATISTICS OF BARKHAUSEN NOISE

In this section, we will address how does our measured value of $\tau = 1.303$ compare with the available data. In nature, numerous systems exhibit crackling noise when driven slowly, characterized by a series of sudden avalanche-like events with a broad range of sizes and distributions. Various studies have been conducted to comprehend the dynamics behind these phenomena, revealing strikingly similar behavior across diverse systems such as earthquakes, plastic deformation, microfractures, vortices in superconductors, dimming events in stars, and Barkhausen noise [8].

Extensive theoretical and experimental investigations, framed within the Langevin theory of domain wall motion, have contributed to the interpretation of Barkhausen noise statistics. These studies have revealed a power law or scaling law distribution aligning seamlessly with experimental data [4, 5].

$$P(s) \sim s^{-\tau} \quad (2)$$

$$P(T) \sim s^{-\alpha} \quad (3)$$

where, s is the jump size and T is the avalanche duration, Durin and Zapperi [6] conducted an investigation into the scaling properties of Barkhausen Noise across several ferromagnetic materials. They measured avalanche distributions, determined scaling constants, and grouped the materials into two classes with τ values of $\tau = 1.50 \pm 0.05$ and $\tau = 1.27 \pm 0.03$. The first category encompasses Si-Fe polycrystals and the

partially crystallized amorphous alloy, while the second category comprises amorphous alloys subjected to stress. These values agree reasonably well with our measured value of τ .

VI. DISCUSSION

Our designed kit serves as a valuable resource for universities and colleges, facilitating the demonstration of the Barkhausen Effect. Students can easily perceive the Barkhausen noise through the integrated speaker. Not only is the kit simple to construct, but it is also highly cost-effective, requiring a modest investment of approximately Rs. 300 (approximately \$3.61 USD). This accessible tool fosters a deeper understanding among students regarding domain walls and the magnetization dynamics of ferromagnetic materials. Additionally, it proves useful for non-destructive testing of materials to identify impurities.

Moreover, the kit enables the calculation of the scaling constant (τ) for the soft iron core, adding an intriguing dimension to its educational utility. The versatility of this kit extends to potential modifications for demonstrating the Barkhausen Effect in various materials. However, it is important to note that the system is susceptible to errors such as mechanical vibrations, mag-

netic field fluctuations, and environmental noise. While data obtained from this device can serve for preliminary calculations, it is acknowledged that more sophisticated and sensitive instruments would be required for detailed analyses. Nevertheless, at the college or school level, this kit proves to be a valuable resource for effective demonstration purposes.

VII. CONCLUSION

In conclusion, the Barkhausen effect is a great educational and research topic. To advance in this field, a foundational understanding is crucial, and this can be achieved through the use of a readily designed kit. Such a kit provides an accessible means for students to delve into the intricacies of the Barkhausen effect, offering valuable insights. The kit is very cost effective and can be introduced in colleges and universities in India and abroad. Moreover, with the aid of this straightforward kit, one can effectively calculate the scaling exponent τ value of a Soft Iron sample and other samples too and can be compared with simulations and other available data.

A. Acknowledgments

I would like to thank *Physical Research Laboratory (PRL) -Ahmadabad* and *Indian Institute of Space Science and Technology (IIST) - Kerala* for providing me all the required facilities.

-
- [1] R. Feynman, Chapter 37: Magnetic materials, *Feynman lectures in physics*, Vol. 2, Publisher: Addison-Wesley (1963).
 - [2] P Weiss , L'hypothèse du champ moléculaire et la propriété ferromagnétique, Vol. 6, Journal: J. Phys. Theor. Appl. (1907).
 - [3] Gianfranco Durin and Stefano Zapperi, The Barkhausen Effect, Vol. 2, Journal: arXiv : [cond-mat/0404512v1](https://arxiv.org/abs/cond-mat/0404512v1) ? (2004).
 - [4] Bruno Alessandro and et. al., Domain-wall dynamics and Barkhausen effect in metallic ferromagnetic materials. I. Theory, Journal: J. Appl. Phys. 68, 2901–2907 (1990).
 - [5] Bruno Alessandro and et. al., Domain-wall dynamics and Barkhausen effect in metallic ferromagnetic materials. II. Experiments, Journal: J. Appl. Phys. 68, 2908–2915 (1990).
 - [6] Gianfranco Durin and Stefano Zapperi, Scaling Exponents for Barkhausen Avalanches in Polycrystalline and Amorphous Ferromagnets, Journal: Physical Review Letters (2000).
 - [7] Sami Kaappa and Lasse Laurson, Barkhausen noise from formation of 360° domain walls in disordered permalloy thin films, Vol. 5, Journal: Physical Review Research (2023).
 - [8] Felipe Bohn and et. al., Playing with universality classes of Barkhausen avalanches, Publisher: Nature- Scientific Reports (2018).
 - [9] Stainslaw and Zurek, Characterisation of Soft Magnetic Materials Under Rotational Magnetisation, Publisher: CRC Press (2019).
 - [10] Irene, LM386 Power Amplifier Circuit, Publisher: ApogeeWeb (2022).